




vista®

for...

COLA'S AND SYRUPS

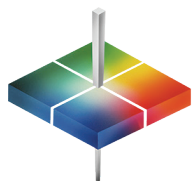

HunterLab

Vista® offers the unique ability to simultaneously measure both visible-range transmittance color and haze% of cola's and syrups, in a single measurement. Vista® includes easy to navigate embedded EasyMatch® Essentials software and modern data communications. With its capacitive high-resolution color touchscreen display, Vista® sets a new standard for transmission color and appearance analysis of cola's and syrups.



Simultaneously measure the transmission color and haze of cola's and syrups



Haze% conforms to ASTM D1003 Procedure B, providing the best correlation to visual assessment



Operator friendly one-touch standardization – no accessories required



Embedded EasyMatch® Essentials software – no PC required




Smart user interface and data management



measure
color
the way your
eye sees it

The Best Solution for Color Measurement of Cola's and Syrups

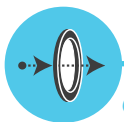

Get instant color and haze results in a single measurement.

Vista® offers this unique ability without the need to change software settings or change the sample container. Both color and haze% can be measured simultaneously in square cuvettes, rectangular cuvettes, and even round vials. This is as simple as it gets.



Smart user interface and data management

View measurement results and color graphics on the 7-inch 1280x800 capacitive industrial touchscreen, and easily communicate data where it needs to go via USB, USB on-the-go, Ethernet, and built-in email.



Embedded EasyMatch® Essentials software – no PC required


The onboard EasyMatch® Essentials software is intuitive and easy to navigate. It offers applications to streamline production, allowing users to create their own workspaces to meet their specific application workflow needs. EasyMatch® Essentials comes preloaded with most standard transmittance scales and indices widely used in industry.

EasyMatch® Essentials




Operator-friendly one-touch standardization

Instrument standardization is fully automated, no tiles or accessories are required. Simply press the 'standardize' icon, and the instrument automatically standardizes for both color and haze.

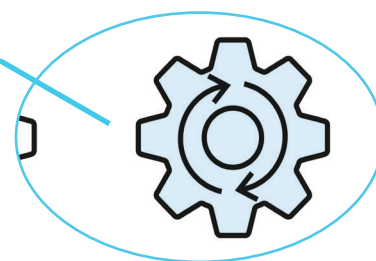

Full spectrum, balanced LED illumination

Illuminating from 400-700nm, the balanced, full spectrum LED light source is guaranteed to perform to original factory specifications for a minimum of 5 years of continuous measurements for years of worry-free operation

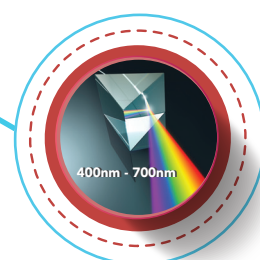

Keys for the successful measurement of Cola's and Syrup



D02-1017-192
Multi-Function Sample Holder



D02-1017-344
Macro Cell Holder for Color
and Haze Measurement

Optional Accessories

L02-1018-683 10mm Flow
Through Cell Holder Kit
Includes:

- 10mm Flow Through Cell (C04-1001-958)
- Cell Holder Assembly (D02-1018-678)
- 1 meter, 7mm ID tubing
- (4) Flow Control Chokes (Various apertures) L02-1018-101

20 mm Flow Through Cell Kit also available

