




vista®

for...

SPIRITS

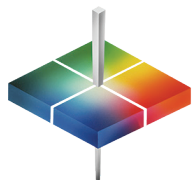

HunterLab

Major brand owners know the importance that color appeal and color consistency play in a consumer's decision to purchase their spirit products. Vista® ensures that the color and clarity of spirits conform to quality standards by simultaneously measuring both transmittance color and haze% in a single measurement. Any issues can be immediately isolated, the root cause discovered, and the process corrected.

Vista® includes easy to navigate embedded EasyMatch® Essentials software and modern data streaming and communications. With its capacitive high-resolution color touchscreen display, Vista® sets a new standard for transmission color and appearance analysis of spirit products.



**Simultaneously measure
transmitted color and haze%
in a single measurement**



**Long-life, full spectrum LED
illumination guaranteed for
five years**




**One-touch standardization -
no standards or accessories
required**



**Embedded EasyMatch® Essentials
software – no PC required**




**Smart user interface and data
management**

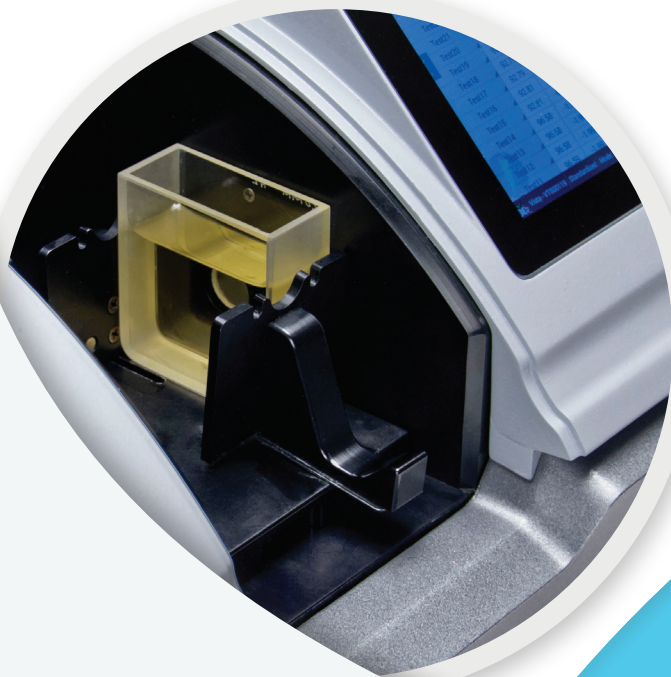


measure
color
the way your
eye sees it

The Best Solution for Color Measurement of Spirits

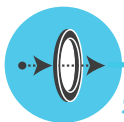

Simultaneously measure transmitted color and haze% in a single measurement

Simultaneously measure transmittance color, and haze% compliant with ASTM D1003 Procedure B, in a single measurement, along with a host of other color scales and color indices used in the beverage industry - in under seven seconds.

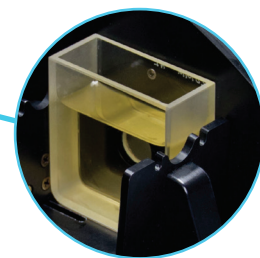

Long-life, full spectrum LED illumination guaranteed for five years

Illuminating from 400-700nm, the balanced, full spectrum LED illumination is guaranteed to perform to original factory specifications for a minimum of 5 years of continuous measurements for years of worry-free.




Operator-friendly one-touch standardization

Instrument standardization is fully automated, no tiles or accessories are required. Simply press the 'standardize' icon, and the instrument automatically standardizes for both color and haze.




Embedded EasyMatch® Essentials software – no PC required

Instrument standardization is fully automated, no tiles or accessories are required. Simply press the 'standardize' icon, and the instrument automatically standardizes for both color and haze.

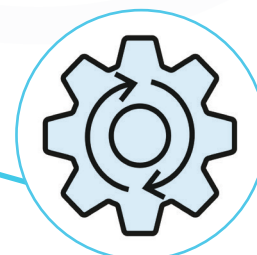

Smart user interface and data management

View measurement results and color graphics on the 7-inch 1280x800 capacitive industrial touchscreen, and easily communicate data where it needs to go via USB, USB on-the-go, Ethernet, and built-in email.

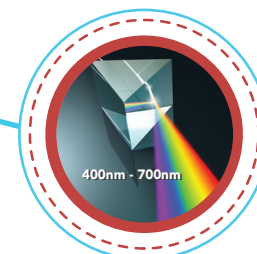

Keys for the successful measurement of Spirits



D02-1017-192
Multi-Function Sample Holder



D02-1017-344
Macro Cell Holder for Color
and Haze Measurement

Optional Accessories

L02-1018-683 10mm Flow Through Cell Holder Kit
Includes:

- 10mm Flow Through Cell (C04-1001-958)
- Cell Holder Assembly (D02-1018-678)
- 1 meter, 7mm ID tubing
- (4) Flow Control Chokes (Various apertures) L02-1018-101

20 mm Flow Through Cell Kit also available

