




vista®

for...

TEA DRINKS

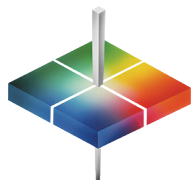

HunterLab

Major brand owners know the importance that color appeal and color consistency play in a consumer's decision to purchase their tea drink products. Vista® ensures that the color and clarity of tea drinks conform to quality standards by simultaneously measuring both transmittance color and haze% in a single measurement. Any issues can be immediately isolated, the root cause discovered, and the process corrected.

Vista® includes easy to navigate embedded EasyMatch Essentials software and modern data streaming and communications. With its capacitive high-resolution color touchscreen display, Vista® sets a new standard for transmission color and appearance analysis of tea drink products.



Simultaneously measure transmitted color and haze% in a single measurement



Long-life, full spectrum led illumination guaranteed for five years




One-touch standardization - no standards or accessories required



Embedded EasyMatch Essentials software – no PC required




Smart user interface and data management

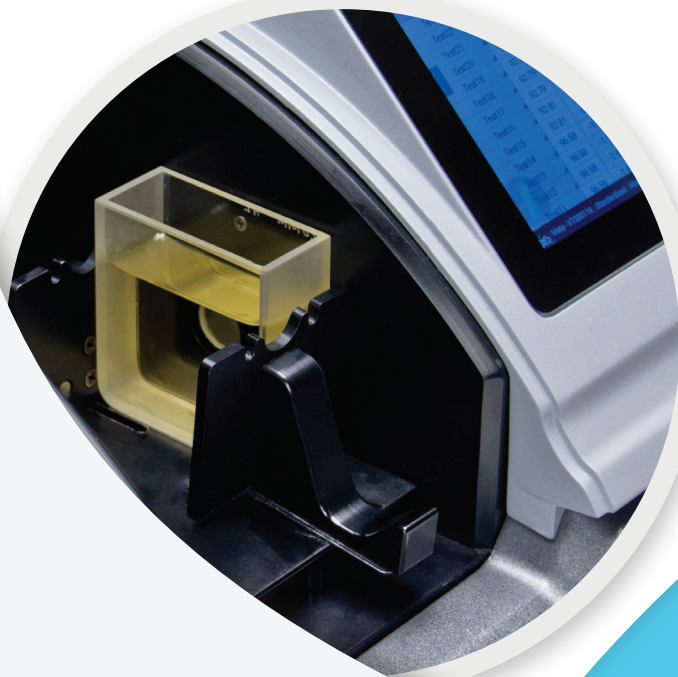


measure
color
the way your
eye sees it

The Best Solution for Color Measurement of Tea Drinks

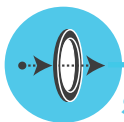

Simultaneously measure transmitted color and haze% in a single measurement

Simultaneously measure transmittance color, and haze% compliant with ASTM D1003 Procedure B, in a single measurement, along with a host of other color scales and color indices used in the beverage industry - in under seven seconds.

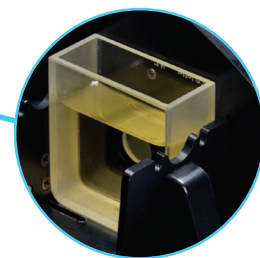

Long-life, full spectrum LED illumination guaranteed for five years

Illuminating from 400-700nm, the balanced, full spectrum LED illumination is guaranteed to perform to original factory specifications for a minimum of 5 years of continuous measurements for years of worry-free.




Operator-friendly one-touch standardization

Instrument standardization is fully automated, no tiles or accessories are required. Simply press the 'standardize' icon, and the instrument automatically standardizes for both color and haze.




Embedded EasyMatch® Essentials software – no PC required

The onboard EasyMatch® Essentials software is intuitive and easy to navigate, offering applications to streamline production and allowing users to create their own workspaces to meet their specific application workflow needs. EasyMatch® Essentials comes preloaded with most standard transmission scales and indices widely used in industry.

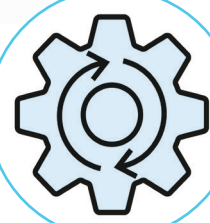

Smart user interface and data management

View measurement results and color graphics on the 7-inch 1280x800 capacitive industrial touchscreen, and easily communicate data where it needs to go via USB, USB on-the-go, Ethernet, and built-in email.

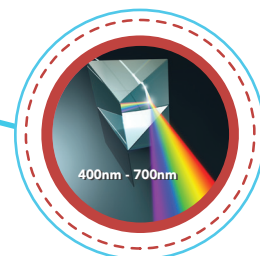

Keys for the successful measurement of Tea Drinks



D02-1017-192
Multi-Function Sample Holder



D02-1017-344
Macro Cell Holder for Color
and Haze Measurement

Optional Accessories

L02-1018-683 10mm Flow Through Cell Holder Kit
Includes:

- 10mm Flow Through Cell (C04-1001-958)
- Cell Holder Assembly (D02-1018-678)
- 1 meter, 7mm ID tubing
- (4) Flow Control Chokes (Various apertures) L02-1018-101

20 mm Flow Through Cell Kit also available

