


vista[®]
Patent Pending

for...

PALM OIL

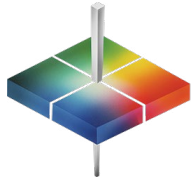

HunterLab

Vista has broken the mold – the 5.25" glass mold that is - when it comes to measuring Palm Oil. Vista includes the unique feature of allowing the user to measure edible oil samples in round vials, or square and rectangular cuvettes, and accurately report standard 1" and 5.25" LOVIBOND[®] RYBN and AOCS RY results. Speed up measurement process time while minimizing product waste and cleanup time. And if product haze analysis is important, Vista simultaneously measures both visible-range transmittance color and haze% in a single measurement.

With modern data communications including USB, Ethernet and Wi-Fi, and a capacitive high-resolution color touchscreen display, Vista sets a new standard for transmission color and appearance analysis of Palm Oil.



Measure in Most Shaped and Sized Vials or Cuvettes, and Accurately Report 1" and 5.25" Lovibond[®] and AOCS Data



Verify Product and Process Integrity - Simultaneously Measure Transmission Color and Haze% in a Single Measurement



Embedded EasyMatch Essentials Quality Control Software - No PC Required



Smart User Interface and Data Management – USB, Ethernet, and Wi-Fi



Operator Friendly One-Touch Standardization – No Accessories Required



measure
color
the way your
eye sees it

The Best Solution for Color Measurement of Palm Oil



Embedded EasyMatch Essentials Quality Control Software – No PC Required

Simultaneously Measure transmittance color and haze% compliant with ASTM D1003 Procedure B, in a single measurement, along with a host of other color scales and color indices in Vista's embedded EasyMatch Essentials quality control software.

Smart User Interface and Data Management

View measurement results and color graphics on the 7-inch 1280x800 capacitive industrial touchscreen, and easily communicate data where it needs to go via USB, USB on-the-go, Ethernet, Wi-Fi, and built in email.




EasyMatch Essentials




Operator-friendly One-touch Standardization

Instrument standardization is fully automated, no tiles or accessories are required. Simply press the 'standardize' icon, and the instrument automatically standardizes for both color and haze.

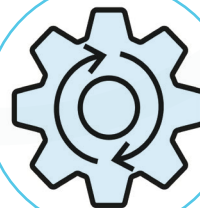

Full Spectrum, Balanced LED Illumination

Illuminating from 400-700 nm, the balanced, full spectrum LED illumination is guaranteed to perform to original factory specifications for a minimum of 5 years of continuous measurements for years of worry-free operation.

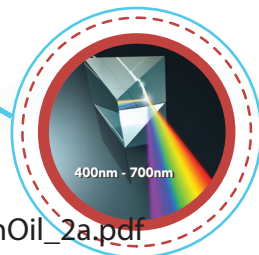

Vista_Solutions_Broch_Food_PalmOil_2a.pdf



Keys for the Successful Measurement of Palm Oil



D02-1017-192
Multi-Function Sample Holder



D02-1017-344
Macro Cell Holder for Color
and Haze Measurement

HunterLab provides a variety of standard cell holder options.
Contact us for help in selecting the best solution to fit your needs.

L02-1018-683
10 mm Flow Through Cell Holder Kit
Includes:

- 10 mm Flow Through Cell (C04-1001-958)
- Cell Holder Assembly (D02-1018-678)
- 1 meter, 7 mm ID tubing
- (4) Flow Control Chokes (Various apertures) (L02-1018-101)

20 mm Flow Through Cell Kit also available

