




vista®

for...

OLIVE OIL

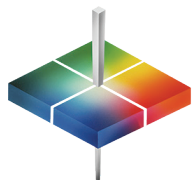


HunterLab

Vista® has broken the mold – the 5.25" glass mold that is - when it comes to measuring olive oil. Vista® includes the unique feature of allowing the user to measure edible oil samples in round vials, or square and rectangular cuvettes, and accurately report standard 1" and 5.25" LOVIBOND® RYBN and AOCS RY results. Speed up measurement process time while minimizing product waste and cleanup time. And if product haze analysis is important, Vista® simultaneously measures both visible-range transmittance color and haze % in a single measurement.

With modern data communications and a capacitive high-resolution color touchscreen display, Vista® sets a new standard for transmission color and appearance analysis of olive oil.



Measure in most shaped and sized vials or cuvettes, and accurately report 1" and 5.25" LOVIBOND® and AOCS data

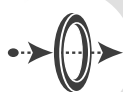

Verify product and process integrity - simultaneously measure transmission color and haze % in a single measurement




Includes easy to navigate embedded EasyMatch® Essentials software – no PC required



Smart user interface and data management




Operator friendly one-touch standardization – no accessories required

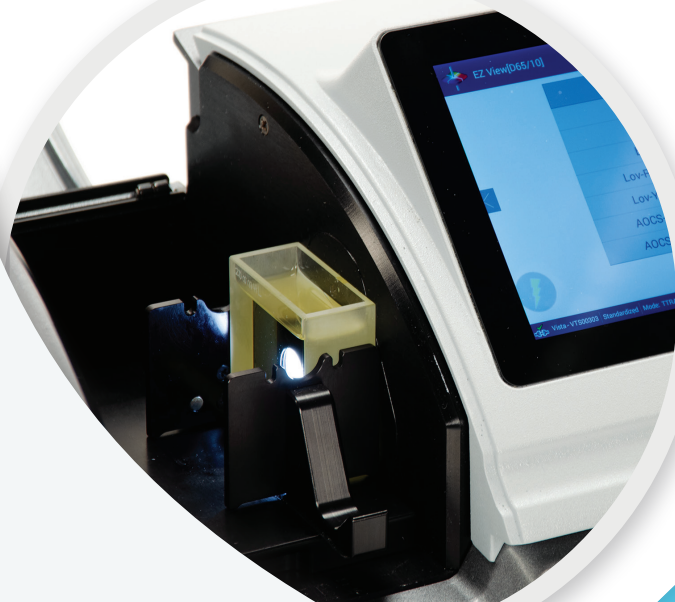

measure
color
the way your
eye sees it

The Best Solution for Color Measurement of Olive Oil

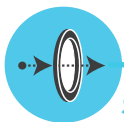

Simultaneously measure both transmittance color, and haze %, in a single measurement

Measure CIE L*, a*, b*, and haze% compliant with ASTM D1003 Procedure B, along with a host of other color scales and color indices in Vista's embedded EasyMatch® Essentials software.



Smart user interface and data management

View measurement results and color graphics on the 7-inch 1280x800 capacitive industrial touchscreen, and easily communicate data where it needs to go via USB, USB on-the-go, Ethernet, and built-in email.




Embedded EasyMatch® Essentials software – no PC required

The onboard EasyMatch® Essentials software is intuitive and easy to navigate. It offers applications to streamline production, allowing users to create their own workspaces to meet their specific application workflow needs. EasyMatch® Essentials comes preloaded with most standard transmittance scales and indices widely used in industry.




Operator-friendly one-touch standardization

Instrument standardization is fully automated, no tiles or accessories are required. Simply press the 'standardize' icon, and the instrument automatically standardizes for both color and haze.

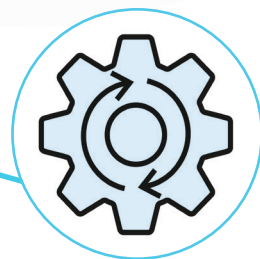

Full spectrum, balanced LED illumination

Illuminating from 400-700nm, the balanced, full spectrum LED illumination is guaranteed to perform to original factory specifications for a minimum of 5 years of continuous measurements for years of worry-free operation.

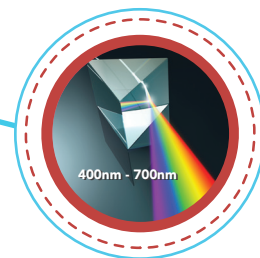

Keys for the successful measurement of Olive Oil



D02-1017-192
Multi-Function Sample Holder



D02-1017-344
Macro Cell Holder for Color
and Haze Measurement

HunterLab provides a variety of standard cell holder options. Contact us for help in selecting the best solution to fit your needs.

Optional Accessories

L02-1018-683 10mm Flow Through Cell Holder Kit Includes:

- 10mm Flow Through Cell (C04-1001-958)
 - Cell Holder Assembly (D02-1018-678)
 - 1 meter, 7mm ID tubing
 - (4) Flow Control Chokes (Various apertures) L02-1018-101
- 20 mm Flow Through Cell Kit also available



LOVIBOND® is a registered trademark of The Tintometer Ltd., UK.